


# Swingaway Lateral Support

# Wheelchairs

This new design of Swingaway Lateral Support enables our standard Lateral Pad to be fitted to a range of wheelchairs with the addition of a swingaway function. The slim design still retains the strength and reliability required to provide consistent, repeatable positioning within a standard wheelchair. They are easy to fit and allow an excellent range of height and width adjustment.

Select the Lateral Pad Size and the Bar Length from the table below. This will indicate the Lateral Offset from the wheelchair backrest tubes. This gives the maximum distance that is possible and will vary as the width between the lateral supports on each side of the wheelchair is increased or decreased.

## Size Information

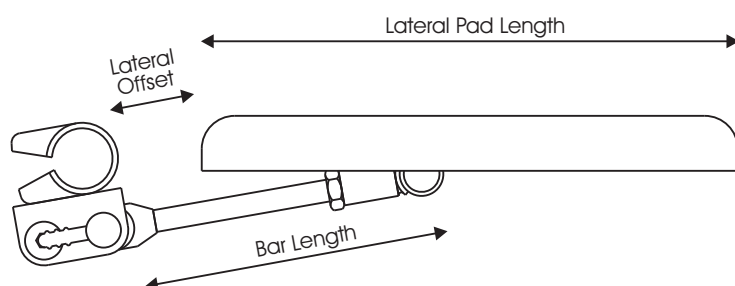

| Part No.   | Part No.   | Lateral Pad Size | LATERAL OFFSET |    |     |      |    |      |    |
|------------|------------|------------------|----------------|----|-----|------|----|------|----|
|            |            |                  | Bar Length     |    |     |      |    |      |    |
| 3/4" Clamp | 7/8" Clamp |                  | 9.5            | 11 | 12  | 13.5 | 16 | 18.5 | 21 |
| C30/21Q    | C30/20Q    | A (13x7.5)       | 1              | 3  | 4   | 5.5  | 8  | 10.5 | 13 |
| C30/21S    | C30/20S    | B (14x7.5)       | 0.5            | 2  | 3   | 4.5  | 7  | 9.5  | 12 |
| C30/21T    | C30/20T    | C (15x9.5)       | 0.5            | 2  | 3   | 4.5  | 7  | 9.5  | 12 |
| C30/21U    | C30/20U    | D (16x9.5)       | 0.5            | 2  | 3   | 4.5  | 7  | 9.5  | 12 |
| C30/21V    | C30/20V    | E (19.5x9.5)     | -              | -  | 0.5 | 2    | 3  | 4.5  | 7  |
| C30/21W    | C30/20W    | F (25x9.5)       | -              | -  | 0.5 | 2    | 3  | 4.5  | 7  |

All dimensions are in cm.

## Compatibility

The Swingaway Lateral Supports are suitable for use with most wheelchairs that have backrest tubes whose diameter is 7/8". They are fitted easily using the split clamp and can be adjusted for width. Height adjustment is provided by moving the lateral support on the wheelchair backrest tubes.

## Additional Support

They can be used to provide a symmetrical trunk position where the support offered by the wheelchair itself is not sufficient.

## Kit Contents

The Swingaway Lateral Support comprises of the following items:

- Lateral support pad
- Pad cover
- 7/8" split tube clamp bracket assembly
- Lateral support bar

## Other Accessories

Lateral Support Pads are also available separately.

They are also available without the swingaway facility

## YOUR PRESCRIPTION DETAILS

LATERAL SUPPORT DETAILS

Lateral Pad Size \_\_\_\_\_

Tube Clamp Size \_\_\_\_\_

Bar Length \_\_\_\_\_

Lateral Offset \_\_\_\_\_

Quantity \_\_\_\_\_

Fabric Colour \_\_\_\_\_

Please attach a copy of this sheet with your order.  
See our Fabric Colour Card for further details of upholstery options.

LFT083-4.cdr

As a company we are committed to promoting a 24 hour approach to postural management through education, design and manufacture. The design of our equipment and training is based on a thorough understanding of the theoretical and practical requirements of the management of posture for people with moderate to severe motor impairment.



**CE** This product is manufactured to comply with the Medical Devices Directive



**Active Design Ltd**  
68K Wyrley Road  
Birmingham B6 7BN  
Tel: (0121) 326 7506  
Fax: (0121) 327 8807  
Web: www.activedesign.co.uk